

Sample preparation – material separation

Urbanti funnels, PMP

Crystal clear. The spiral-shaped ribs increase the rate of filtration and prevent the trapping of air between the filter paper and the funnel. With long stem.


Volume approx. ml	Ø mm	Length mm	Stem Ø mm	Stem length mm	PU	Cat. No.
30	51	195	3	150	6	325095
80	70	210	3	150	6	325195
250	100	198	7	108	4	325295
630	140	247	10	132	3	325395
1800	196	315	20	155	2	325495



Analytical funnels, PP

Transparent. With long stem and grooves. Rigidified by a thickened edge. Rapid flow due to a steep 60° angle.

Volume approx. ml	Ø mm	Length mm	Stem Ø mm	Stem length mm	PU	Cat. No.
50	50	194	5	150	10	80162
100	72	208	5	143	10	80164
225	91	227	5	145	10	80165

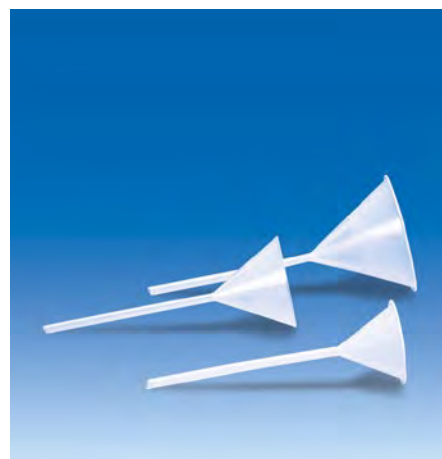

Büchner funnels, PP

Three parts. Upper and lower parts are detachable to facilitate cleaning.

Volume approx. ml	Filter Ø mm	Length mm	Hole Ø mm	PU	Cat. No.
40	42.5	95	1.2	1	80437
70	55	113	1.1	1	80438
180	70	145	2.0	1	80439
280	80	165	2.0	1	80440
390	90	180	2.5	1	80441
810	110	210	2.5	1	80442
2100	160	280	2.75	1	80443
6000	240	350	3.0	1	80445

